


# REVERSE BEAEBSE PEEL



Direct line to contact dermatologists and research biologists  
for scientific questions or enquiries regarding outpatient treatment:  
[scientificdiv.gpq@wiqo.eu](mailto:scientificdiv.gpq@wiqo.eu)

GPQ srl

Operational Headquarters: Via Flavia, 23/1 - 34148 Trieste - Italy  
Tel.: +39 040 360408 - Fax: +39 040 761393 - [gpq@wiqo.eu](mailto:gpq@wiqo.eu)  
[www.wiqo.com](http://www.wiqo.com)

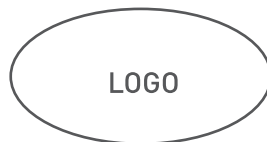

DISTRIBUTOR  
contacts  
contacts

## PROTOCOL FOR USE



## INDICATIONS

Treatment of **melasma** and **skin hyperpigmentation**.

## FREQUENCY OF USE

Every 7 to 30 days, according to the doctor's assessment. Start with at least 3 weekly sessions and then reduce the frequency, according to the intensity of melasma and its improvement.

## METHOD OF USE AND APPLICATION



RsP 1

- Clean the area to be treated with WiQo **P SOLUTION**.
- Apply generously to all affected areas and massage delicately with nitrile gloves until absorbed.
- Repeat the application a second time.
- Once completely absorbed, wash the area with a soaked pad.
- Dry well.



RsP 2

- Apply superficially with the nitrile gloves on pigmented areas only, without massaging it in.
- Leave on for 3 minutes.\*
- If the patient reports a burning sensation higher than 5/10 or if reddening occurs, apply the **NEUTRALISER** immediately.



N

- Shake the **NEUTRALISER** before use.
- Apply the product. Neutralisation takes place over a few seconds.
- If the patient reported burning, leave to rest until burning decreases.
- Remove with a pad soaked with water.
- Dry well.



RsP 3

- Apply superficially with the nitrile gloves on pigmented areas only, without massaging it in.
- Leave on for 3 minutes.\*
- If the patient reports a burning sensation higher than 5/10 or if reddening occurs, apply the **NEUTRALISER** immediately.



N

- Shake the **NEUTRALISER** before use.
- Apply the product. Neutralisation takes place over a few seconds.
- If the patient reported burning, leave to rest until burning decreases.
- Remove with a pad soaked with water.
- Dry well.



## IMMEDIATE POST-TREATMENT AND HOME TREATMENT

- Apply WiQo **REVERSE SERUM** with a delicate massage until its absorption.
- Apply WiQo **ICP CREAM** or another sun protection cream (SPF 50+).

## CONTRAINDICATIONS

- Irritated skin
- Irritable skin
- Seborrhoeic dermatitis
- Atopic dermatitis
- Pregnancy and breastfeeding
- Active skin rashes and/or herpetic lesions

## WARNINGS

- For medical professional use only.
- Medical device for external use.
- **The application of Phase 2 and Phase 3 must be performed with care:**
  - If redness or strong burning sensation occurs during application, immediately apply Neutraliser.
  - Do not apply Phase 3 if the skin is highly reddened or the patient reports a persistent burning sensation after the neutralisation of Phase 2.
  - If raised spots or small off-white areas appear on treated areas, neutralise immediately.
  - If the burning sensation or an intense rash persists after using the Neutraliser and subsequently rinsing the skin, apply a steroid cream.
- Do not use latex gloves as they may cause allergies. Nitrile gloves are recommended.
- Treat following a reasonable time after shaving or waxing.
- Avoid contact with eyes. In case of contact with eyes rinse thoroughly with saline solution, to be kept at hand at all time.
- Avoid contact with mucous membranes. In case of contact with mucous membranes, immediately apply the Neutraliser and then remove with a soaked pad.
- The product is corrosive.

\* RsP 2 and RsP 3 can be left on for up to 15 minutes in cases of very thick skin.